

*Ostomy Care Product Catalogue*

# ***Eliminating*** the burden of leakage



**Coloplast**



*“Ideally, we would like to have people forget they have a medical condition. It’s all about people just living the life they want.”*

Kristian Villumsen  
CEO Coloplast

## SenSura® Mio

- 4 SenSura® Mio  
SenSura® Mio Convex

## SenSura® Mio Baby and Kids

- 55 Long wear time - less disturbance so kids can move like kids

## SenSura®

- 20 **Dual protection - less leakage, healthier skin**  
Double Protection Adhesive  
SenSura® Closed Bag  
SenSura® Drainable Bag  
SenSura® Urostomy Bag  
SenSura® Post-Op Bag  
SenSura® Baseplates

## Alterna® Free

- 28 **Spiral security with more freedom**  
Alterna® Free Spiral Adhesive  
Alterna® Free Closed Bag  
Alterna® Free Drainable Bag  
Alterna® Free Urostomy Bag  
Alterna® Free Baseplates

## Hospital Assortment

- 35 SenSura® Mio  
SenSura®  
Alterna®

## Brava®

- 43 **For customising appliance fit and for issues such as leakage and skin irritation.**  
Brava® Adhesive Remover Spray and Wipes  
Brava® Powder  
Brava® Barrier Cream  
Brava® Belt  
Brava® Mouldable Ring  
Brava® Paste (Tube)  
Brava® Strip Paste  
Brava® Elastic Tape  
Brava® Protective Sheet  
Brava® Skin Barrier Spray and Wipes  
Brava® Lubricating Deodorant  
Brava® Protective Seal  
Brava® Belt for SenSura® Mio

## Fistula and Wound Management System

- 57 Composition



SenSura<sup>®</sup> Mio

***Follows shapes  
and movement***

The fit is the ***difference***




For more ***than 60 years***, Coloplast has been developing products and services to ***make life easier*** for people with personal and private medical conditions.


SenSura® Mio ostomy appliance with its elastic adhesive is a milestone in our mission.



# Leakage is the *number 1* problem

Leakage is *the number one problem* for all people with a stoma.<sup>1</sup>

91%   
worry about  
leakage <sup>2</sup>


80%   
say they have  
leakage issues <sup>2</sup>

Poor appliance fit to body is the main cause of leakage and skin irritation. Getting a secure fit is therefore key to preventing leakage and skin issues.

## *The right fit* makes all the difference

Finding the right fit is an important part of improving stoma-related quality of life. However, body shapes, stomas and the location of the stoma opening vary greatly from person to person. These factors can make it hard to find the right fit for a patient's body profile.

Users can be divided into three distinct body profiles.

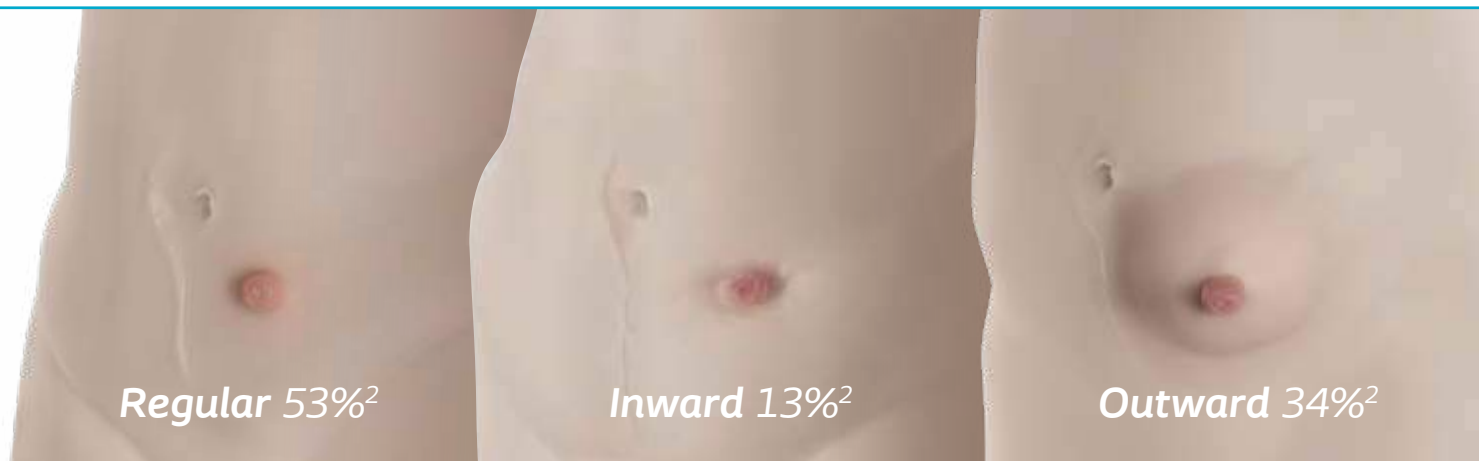

# Introducing SenSura® Mio with BodyFit Technology®



BodyFit  
Technology®



A common fit challenge for users with regular body profiles is securing a good fit that supports users' lifestyles, even physically active ones with their associated abdominal movement.<sup>3</sup>

SenSura® Mio's unique elastic adhesive follows body movements. Whether the patient is bending down, stretching or twisting, this appliance maintains a secure fit.

The baseplate is designed with BodyFit Technology that secures the right fit for your patient's body, no matter their body shape and specific challenges.

Follows  
shapes and  
movements



**86%**  
experience a  
good or very  
good fit<sup>4</sup>






# The science behind SenSura<sup>®</sup> Mio *elastic adhesive*

A special blend of hydrocolloids, polymers and tackifiers are the main ingredients of the SenSura Mio elastic adhesive formula.

1 Special blend of hydrocolloids, polymers and tackifiers

2 Elastic topfilm

 Hydrocolloid  
 Polymer

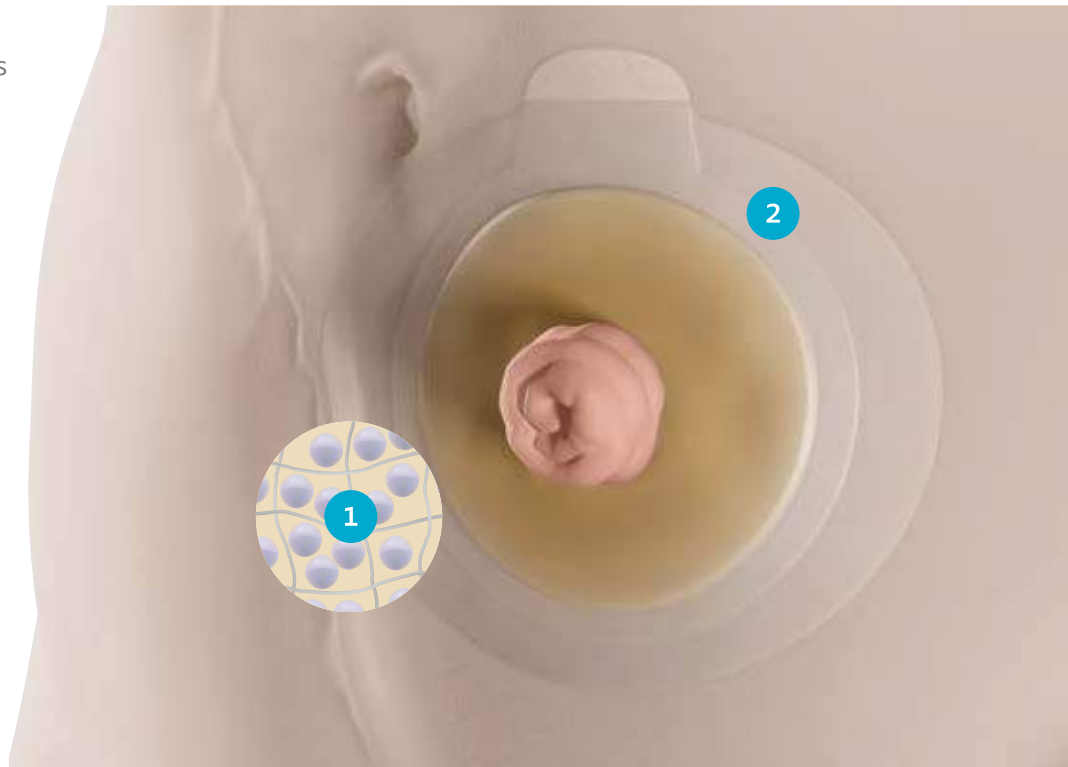

Thanks to **BodyFit Technology<sup>®</sup>**, the baseplate can bend and stretch with the body while maintaining a secure fit to the skin.



**2x**  
more stretchable  
than a conventional  
appliance<sup>5</sup>

# 3 more reasons why SenSura<sup>®</sup> Mio is *a great choice for your patients*



## Discreet bag colour

The all-textile material in specially chosen, neutral colours is designed to be unnoticeable, even under white clothing.

## A filter you can trust

A unique, full-circle filter provides greater protection against ballooning.

Proven to reduce incidents of  
ballooning by **61%**<sup>6</sup>




## A choice of couplings

**Click coupling** has a lock ring designed for reassurance and a feeling of security. An audible 'click' confirms when the lock ring is correctly locking the bag to the baseplate.

*SenSura<sup>®</sup> Mio Convex*

# *Convex without compromise*

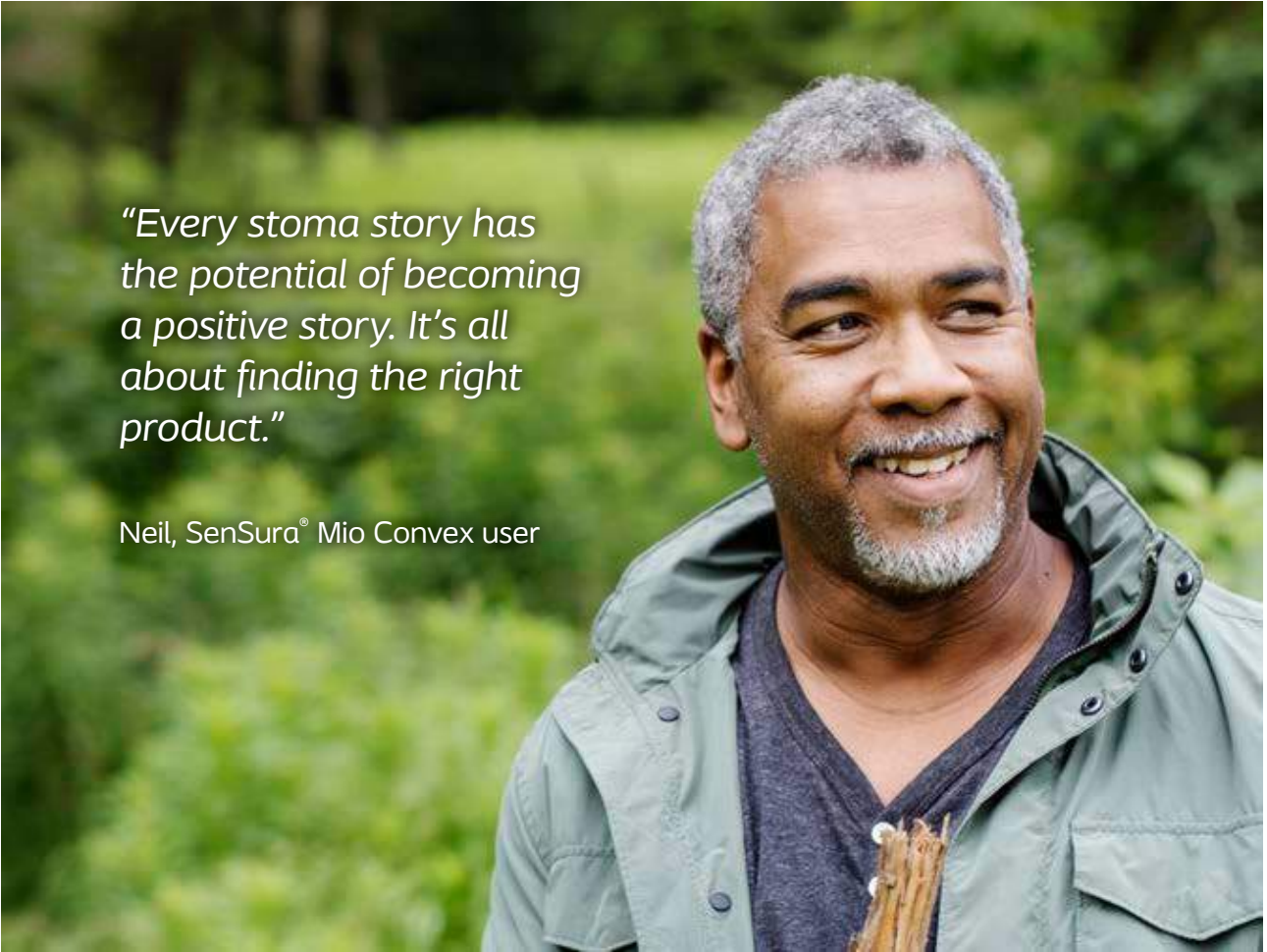

A photograph of a middle-aged Black man with grey hair and a goatee, smiling and looking slightly to his right. He is wearing a green jacket over a dark shirt. The background is a blurred green field.

*"Every stoma story has the potential of becoming a positive story. It's all about finding the right product."*

Neil, SenSura® Mio Convex user

For more ***than 60 years***, Coloplast has been developing products and services to ***make life easier*** for people with personal and private medical conditions.

SenSura® Mio ostomy appliance with its elastic adhesive is a milestone in our mission.

# Introducing SenSura<sup>®</sup> Mio with BodyFit Technology<sup>®</sup>



Proven  
leakage  
reduction<sup>1</sup>

Some fit challenges for users with inward body profiles include:

- Stomas are commonly level with or below skin surface<sup>2</sup>
- High risk of leakage (84%)<sup>2</sup>
- Traditional convex appliances can feel bulky and stiff<sup>3</sup>

For too long, it has been accepted that finding a more secure convex appliance means settling for less flexibility.

SenSura<sup>®</sup> Mio Convex provides flexibility for people with stomas that are retracted, flush or in an inward area – without compromising security.



Traditional appliances  
Potential leakage

SenSura<sup>®</sup> Mio Convex comes in three variants: **soft**, **light** and **deep**. The deep Convex is twice as flexible as a traditional light Convex.<sup>4</sup>



Soft



Light



Deep

## 86%

Assessed the level of comfort when wearing the baseplate as good or very good<sup>5</sup>



# Proven security for inward areas

The SenSura® Mio Convex adhesive features BodyFit Technology – an adhesive technology that secures the right fit for every body, no matter their body shape and challenges.

The adhesive's patented, integrated flexlines give the flexibility needed for comfortable bending, without compromising on stability.

The material in the flexlines is slightly thinner than the surrounding material, and their pattern ensures the convex bends around the stoma – not in the middle of it.



A  
comfortable  
convex









## 2x

Even the deep convex is twice as flexible as the traditional light convex<sup>4</sup>












































# Sensura<sup>®</sup> Mio

## *the fit is the difference*

One-piece	<p>Sensura<sup>®</sup> Mio</p> 	<p>Sensura<sup>®</sup> Mio Convex</p>  
Two-piece, Click	<p>Sensura<sup>®</sup> Mio</p> 	<p>Sensura<sup>®</sup> Mio Convex</p>  

References: 1. Coloplast, Market research, OLS 2016\_ECET Oral Presentation, Data on file PM-05164. 2. Claessens et. al. Gastrointestinal Nursing 2015. 3. Coloplast Market study, Ostomy Life Study 2016, data on file PM-05235. 4. Coloplast, Product Evaluation, SenSura Mio, 2014 (VV-0197319). 5. Coloplast, Laboratory Report, Stretchability of SenSura Mio (VV-0197750). 6. Coloplast, Clinical Study, 2010 (VV-0009690) for ileostomists, compared to Standard of Care

Product Code	Product Description	Pack Size
<b>SenSura® Mio One Piece Closed Bag</b>		
108710	SenSura® Mio One Piece, Soft Flat Base Plate, Closed, Midi 15-45mm, Opaque with Inspection Window	30
108810	SenSura® Mio One Piece, Soft Flat Base Plate, Closed, Maxi 15-45mm, Opaque with Inspection Window	30
108910	SenSura® Mio One Piece, Light Flat Base Plate, Closed, Midi 15-55mm, Opaque with Inspection Window	30
163010	Sensura® Mio One Piece, Soft Convex Base Plate, Closed, Midi 15-33mm, Opaque with Inspection Window	10
163060	Sensura® Mio One Piece, Soft Convex Base Plate, Closed, Maxi 10-50mm, Opaque with Inspection Window	10
163220	Sensura® Mio One Piece, Light Convex Base Plate, Closed, Midi 15-33mm, Opaque with Inspection Window	10
163230	SenSura® Mio One Piece, Light Convex Base Plate, Closed, Midi 31-43mm, Opaque with Inspection Window	10
163320	SenSura® Mio One Piece, Light Convex Base Plate, Closed, Maxi 10-43mm, Opaque with Inspection Window	10
163420	SenSura® Mio One Piece, Deep Convex Base Plate, Closed, Midi 15-33mm, Opaque with Inspection Window	10
163430	SenSura® Mio One Piece, Deep Convex Base Plate, Closed, Midi 31-43mm, Opaque with Inspection Window	10
163520	SenSura® Mio One Piece, Deep Convex Base Plate, Closed, Maxi 10-43mm, Opaque with Inspection Window	10
<b>SenSura® Mio One Piece Drainable Bag</b>		
104310	SenSura® Mio One Piece, Flat Base Plate, Drainable, Mini , Opaque with Inspection Window, 10-45mm	30
104410	SenSura® Mio One Piece, Flat Base Plate, Drainable, Midi Opaque with Inspection Window, 10-45mm	30
104510	SenSura® Mio One Piece, Flat Base Plate, Drainable, Maxi Opaque with Inspection Window, 10-55mm	30
104560	SenSura® Mio One Piece, Flat Base Plate, Drainable, Maxi Opaque , 10-55mm	30
104580	SenSura® Mio One Piece, Flat Base Plate, Drainable, Maxi Transparent, 10-55mm	30
164010	SenSura® Mio One Piece, Soft Convex Base Plate, Drainable Midi 15-33mm, Opaque with Inspection Window	10
164060	SenSura® Mio One Piece, Soft Convex Base Plate, Drainable, Maxi 10-50mm, Opaque with Inspection Window	10
164150	SenSura® Mio One Piece, Soft Convex Base Plate, Drainable, Maxi Transparent 10-50mm	10
164210	SenSura® Mio One Piece, Light Convex Base Plate, Drainable, Midi 10-23mm, Opaque with Inspection Window	10
164220	SenSura® Mio One Piece, Light Convex Base Plate, Drainable, Midi 15-33mm, Opaque with Inspection Window	10
164230	SenSura® Mio One Piece, Light Convex Base Plate, Drainable, Midi 31-43mm, Opaque with Inspection Window	10
164350	SenSura® Mio One Piece, Light Convex Base Plate, Drainable, Maxi Transparent 10-43mm	10
164360	SenSura® Mio One Piece, Light Convex Base Plate, Drainable, Maxi 10-43mm, Opaque with Inspection Window	10
164410	SenSura® Mio One Piece, Deep Convex Base Plate, Drainable, Midi 10-23mm, Opaque with Inspection Window	10
164420	SenSura® Mio One Piece, Deep Convex Base Plate, Drainable, Midi 15-33mm, Opaque with Inspection Window	10
164430	SenSura® Mio One Piece, Deep Convex Base Plate, Drainable, Midi 31-43mm, Opaque with Inspection Window	10
164550	SenSura® Mio One Piece, Deep Convex Base Plate, Drainable, Maxi Transparent 10-43mm	10
164560	SenSura® Mio One Piece, Deep Convex Base Plate, Draianble, Maxi 10-43mm, Opaque with Inspection Window	10
<b>SenSura® Mio One Piece Urostomy</b>		
105810	SenSura® Mio One Piece, Flat Base Plate, Urostomy, Midi Opaque , 10-45mm	30
105820	SenSura® Mio One Piece, Flat Base Plate, Urostomy, Maxi Opaque, 10-45mm	30
105860	SenSura® Mio One Piece, Flat Base Plate, Urostomy, Maxi Transparent, 10-45mm	30
168010	SenSura® Mio One Piece, Soft Convex Base Plate, Urostomy, Midi Opaque 10-33mm	10
168050	SenSura® Mio One Piece, Soft Convex Base Plate, Urostomy, Maxi Opaque 10-33mm	10
168100	SenSura® Mio One Piece, Soft Convex Base Plate, Urostomy, Maxi Transparent 10-50mm	10
168220	SenSura® Mio One Piece, Light Convex Base Plate, Urostomy, Midi Opaque 15-33mm	10
168230	SenSura® Mio One Piece, Light Convex Base Plate, Urostomy, Midi Opaque 31-43mm	10
168260	SenSura® Mio One Piece, Light Convex Base Plate, Urostomy, Maxi Opaque 15-33mm	10
168270	SenSura® Mio One Piece, Light Convex Base Plate, Urostomy, Maxi Opaque 31-43mm	10
168370	SenSura® Mio One Piece, Light Convex Base Plate, Urostomy,Maxi Transparent 10-43mm	10
168520	SenSura® Mio One Piece, Deep Convex Base Plate, Urostomy, Midi Opaque 15-33mm	10
168570	SenSura® Mio One Piece, Deep ConvexBase Plate, Urostomy, Maxi Opaque 31-43mm	10
168670	SenSura® Mio One Piece, Deep Convex Base Plate, Urostomy, Maxi Transparent 10-43mm	10
<b>SenSura® Mio CLICK Two Piece Base Flat</b>		
105020	SenSura® Mio Flat Base Plate, Two piece, Click, With Belt Ears, 40mm, 10-35	5 
105120	SenSura® Mio Flat Base Plate, Two piece, Click, With Belt Ears, 50mm, 10-45	5 
105220	SenSura® Mio Flat Base Plate, Two piece, Click, With Belt Ears, 60mm, 10-55	5 
105310	SenSura® Mio Flat Base Plate, Two piece, Click, With Belt Ears, 70mm, 10-65	5 
<b>SenSura® Mio CLICK Two Piece Base Plate Convex</b>		
169010	SenSura® Mio Two Piece, Click 40mm Convex Light 10-20 base plate	5 
169110	SenSura® Mio Two Piece, Click 50mm Convex Light 15-30 base plate	5 
169210	SenSura® Mio Two Piece, Click 60mm Convex Light 15-40 base plate	5 
169310	SenSura® Mio Two Piece, Click 70mm Convex Light 15-53 base plate	5 
169410	SenSura® Mio Two Piece, Click 40mm Convex Deep 10-20 base plate	5 
169510	SenSura® Mio Two Piece, Click 50mm Convex Deep 15-30 base plate	5 
169610	SenSura® Mio Two Piece, Click 60mm Convex Deep 15-40 base plate	5 
169710	SenSura® Mio Two Piece, Click 70mm Convex Deep 15-53 base plate	5 

Product Code	Product Description	Pack Size
<b>SenSura® Mio Two Piece Closed Bag</b>		
114020	SenSura® Mio Two Piece, Closed, Click, Maxi Opaque with Inspection-Window, 40mm	30 
114030	SenSura® Mio Two Piece, Closed, Click, Maxi Transparent, 40mm	30 
114130	SenSura® Mio Two Piece, Closed, Click, Maxi Opaque with Inspection-Window, 50mm	30 
114140	SenSura® Mio Two Piece, Closed, Click, Maxi Transparent, 50mm	30 
114220	SenSura® Mio Two Piece, Closed, Click, Maxi Opaque with Inspection-Window, 60mm	30 
114230	SenSura® Mio Two Piece, Closed, Click, Maxi Transparent, 60mm	30 
114250	SenSura® Mio Two Piece, Closed, Click, Maxi Opaque with Inspection-Window, 70mm	30 
114260	SenSura® Mio Two Piece, Closed, Click, Maxi Transparent, 70mm	30 
<b>SenSura® Mio CLICK Two Piece Drainable Bag</b>		
114310	SenSura® Mio, Two-piece, Click, Drainable, Midi Window, 40mm	30 
114320	SenSura® Mio, Two-piece, Click, Drainable, Maxi Window, 40mm	30 
114330	SenSura® Mio, Two-piece, Click, Drainable, Maxi Transparent, 40mm	30 
114340	SenSura® Mio, Two-piece, Click, Drainable, Maxi Opaque, 40mm	30 
114420	SenSura® Mio, Two-piece, Click, Drainable, Midi Window, 50mm	30 
114430	SenSura® Mio, Two-piece, Click, Drainable, Maxi Window, 50mm	30 
114440	SenSura® Mio, Two-piece, Click, Drainable, Maxi Transparent, 50mm	30 
114450	SenSura® Mio, Two-piece, Click, Drainable, Maxi Opaque, 50mm	30 
114510	SenSura® Mio, Two-piece, Click, Drainable, Midi Window, 60mm	30 
114520	SenSura® Mio, Two-piece, Click, Drainable, Maxi Window, 60mm	30 
114530	SenSura® Mio, Two-piece, Click, Drainable, Maxi Transparent, 60mm	30 
114540	SenSura® Mio, Two-piece, Click, Drainable, Maxi Opaque, 60mm	30 
114560	SenSura® Mio, Two-piece, Click, Drainable, Maxi Window, 70mm	30 
114570	SenSura® Mio, Two-piece, Click, Drainable, Maxi Transparent, 70mm	30 
<b>SenSura® Mio CLICK Two Piece Urostomy Bag</b>		
114910	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Opaque 40mm	30 
114920	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Transparent 40mm	30 
114960	SenSura® Mio, Two-piece, Click, Urostomy, Midi Opaque 50mm	30 
114970	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Opaque 50mm	30 
114980	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Transparent 50mm	30 
114990	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Opaque 60mm	30 
115000	SenSura® Mio, Two-piece, Click, Urostomy, Maxi Transparent 60mm	30 
<b>SenSura® Mio Kids</b>		
187110	SenSura® Mio Kids, One Piece, Midi Drainable bag, Transparent, Cut to fit 0-35mm	10
187120	SenSura® Mio Kids One Piece Urostomy bag, Transparent, Cut to fit 0-35mm	10
187210	SenSura® Mio Kids Two Piece, Mini , Drainable Bag, Transparent 27mm	10
187230	SenSura® Mio Kids Two Piece, Midi , Drainable Bag, Transparent 27 mm	10
187250	SenSura® Mio Kids Two Piece, Midi Drainable Bag, Transparent 35mm	10
187300	SenSura® Mio Kids Two Piece, Urostomy Bag, Opaque 27mm	10
187310	SenSura® Mio Kids Two Piece, Urostomy Bag, Opaque 35mm	10
187350	SenSura® Mio Kids Two Piece, Base Plate CM27mm	5
187470	SenSura® Mio Kids Two Piece, Base Plate CM35mm	5
<b>SenSura® Mio Baby</b>		
187000	SenSura® Mio Baby Two Piece, Sterile Drainable Bag, Transparent	10
187050	SenSura® Mio Baby Two Piece Baseplate 0-40mm	5
120140	Brava Baby Adhesive remover	30

# Introducing Coloplast® Baby

*Long wear time – less disturbance*

**Coloplast® Baby**

**Coloplast Baby**  
delivers long wear time  
without compromising  
skin-friendliness.



Long  
wear time<sup>1</sup>

# The *first* appliance designed specially for premature babies and infants

Coloplast® Baby is a two-piece appliance. It features an elastic adhesive that balances the pH level of the skin and a sterile ostomy bag. Crucially, it provides a secure fit to body and is skin-friendly.



## Coloplast Baby – Benefits for professionals

Coloplast Baby is designed to ease your daily practice and:

- Reduce the time spent handling the infant
- Minimise disturbance for the infant
- Makes it easier to meet individual needs (no starter hole)
- Enable access to the stoma for refeeding

## Additional infant-friendly features



### Elastic adhesive

Coloplast Baby is designed to fit individual body profiles and cover multiple stomas. Brava® Baby Adhesive Remover shall be used when changing the adhesive.



### Soft flange

The coupling system includes a soft flange so the bag can be attached to the baseplate with little pressure on the abdomen.



### Soft, easy-to-use outlet

The outlet is compatible with enteral syringes for easy emptying while the soft material is gentle on baby skin.



### Refeeding window

The first dedicated refeeding window on an ostomy appliance for easy facilitation of refeeding and access to the stoma, for less disturbance of the baby.

## Product assortment list

Product name	Item code
Coloplast Baby 2p Flex bag	18700
Coloplast Baby 2p Flex baseplate	18705
Brava Baby Adhesive Remover	12014

# Introducing Coloplast® Kids

*So kids can move like kids*

**Coloplast® Kids**

**Coloplast® Kids**  
is designed to fit the  
contours of kids' small,  
round bellies.






# The first appliance designed specially for kids

Coloplast® Kids features an elastic adhesive, enabling it to adapt to individual body movements. It's available in both a 1-piece and 2-piece solution.

## Additional kid-friendly features

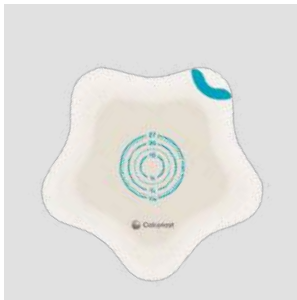


### Star-shaped elastic adhesive

The adhesive is skin-friendly and star-shaped with special petals designed to hug kids' small, round bellies while giving more flexibility in terms of placement on the abdomen.



### Fit zones

Special fit zones in the rim and centre of the appliance increase its flexibility around the stoma for optimal range of movement.

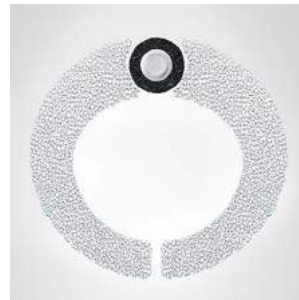

### Unique, full-circle pre-filter

This specially developed feature is proven to postpone clogging. It also reduces incidents of ballooning by 61%.<sup>1</sup>



### Discreet material

Coloplast® Kids features a neutral-coloured, water-repellent, comfortable textile for a non-medical, discreet appearance.

## Product assortment list

	Item code	Product name
Kids 1P	18711	Coloplast Kids 1p op md trp 0-35
	18712	Coloplast Kids 1p uro 10-35 trp
Kids 2P bags	18721	Coloplast Kids 2p Flex op mn 27mm trp
	18723	Coloplast Kids 2p Flex op md 27mm trp
	18725	Coloplast Kids 2p Flex op md 35mm trp
	18730	Coloplast Kids 2p Flex uro 27mm
	18731	Coloplast Kids 2p Flex uro 35mm
Kids 2P baseplates	18735	Coloplast Kids 2p Flex 27mm, 0-27
	18747	Coloplast Kids 2p Flex 35mm, 0-35

## Want to learn more?

Learn more about babies and kids with "The Global Paediatric Stoma Care Best Practice Guideline" developed by our Global Paediatric Stoma Care Nurses Advisory Board. Get a free copy at [Coloplast® professional.com](https://www.coloplastprofessional.com)

Coloplast®  
**Professional**

1. Coloplast, Clinical Study, 2010, 61% for ileostomists, 53% for colostomists.



# *SenSura<sup>®</sup>*

*Dual protection  
- less leakage,  
healthier skin*

**SenSura<sup>®</sup>**


# Double-protection adhesive

Ostomates need an appliance that provides optimal security and skin protection. With SenSura® they have a solution they can trust.

The unique SenSura® double-layer adhesive protects the peristomal skin throughout the entire wear time.

## Protection layer


- Strong erosion-resistant adhesive protects against stoma effluent
- Creates a tight seal between adhesive and stoma

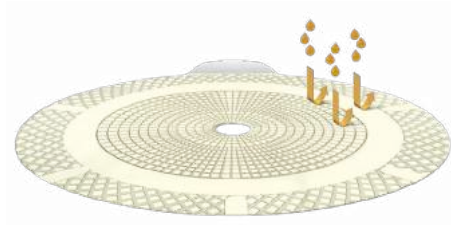

Protects against stoma output

## Skin care layer

- Absorbs excess moisture to keep peristomal skin healthy
- Ensures consistent adhesion

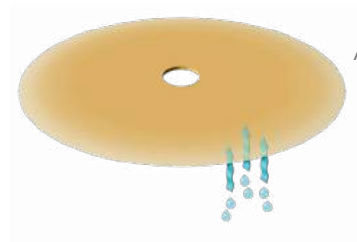

Absorbs excess moisture



**Together the two layers of the double-protection adhesive offer security and skin care**

---

## Two variants of the double-protection adhesive

### SenSura®

The original adhesive for optimal skin friendliness.



### SenSura® Xpro

For extra protection against aggressive output.



**The combination of nursing care and SenSura® double-layer adhesive is proven to:**

- **Reduce** incidents of leakage
- **Improve** peristomal skin conditions
- **Improve** quality of life

# SenSura® Closed Bags

## SenSura® 1-Piece Closed

SenSura® 1-piece closed has an integrated filter which avoids the risk of the pouch ballooning up and causing unnecessary bulges, and neutralises odour. The removal ear makes it easy to remove the adhesive without leaving residue on the skin, allowing for more comfortable pouch changes. The soft backing fabric of SenSura® Closed is strong and water repellent, allowing for easy drying after a wash or swim.



## SenSura® 2-Piece Closed

### SenSura® Flex

SenSura® Flex Closed is fixed securely to the SenSura double-layer adhesive baseplate by a strong adhesive coupling system. The adhesive coupling is flexible, allowing freedom of movement and making it comfortable and discreet to wear. The pouch can be removed from the SenSura® double-layer baseplate, allowing the pouch to be changed more often than the adhesive baseplate.



### SenSura® Click

SenSura® Click 2-piece combines the unique double-layer adhesive with a removable ostomy bag. The bag is connected to the base plate by a mechanical click-lock coupling system which is leak proof and secure. The simplicity of SenSura® Click 2-Piece enables minimal pressure and can be removed with an easy “up-and-over” motion.



# SenSura® Drainable Bags

## ***SenSura® 1-Piece Drainable with Hideaway***

To change the pouch, the whole appliance is removed and replaced. SenSura® 1-piece drainable is available with flat or convex baseplates. The broad assortment consists of various pre-cut or hole sizes which can be customised, with pouches in a variety of sizes available in transparent or opaque colour.



## ***SenSura® 2-Piece Drainable with Hideaway***

### ***SenSura® Flex***

SenSura® Flex has a soft, flexible coupling combined with a floating flange that adjusts to individual body contours, while ensuring a stable position of the pouch. It has the flexibility of a one-piece, while providing the advantages of a two-piece. The secure coupling system is conveniently flat without hard elements, making it discreet and easy to handle.



### ***SenSura® Click***

A secure two-piece mechanical coupling, with an audible click for extra reassurance. SenSura® Click provides a gentle application by minimising pressure on the abdomen, through smooth coupling of pouch to the baseplate. A click sound provides an audible guarantee for secure locking of the pouch to the baseplate.





# SenSura® Urostomy Bags

## **SenSura® 1-Piece Urostomy**

With SenSura® 1-piece urostomy, the SenSura® double-layer adhesive is permanently fixed to the pouch. To change the pouch, the whole appliance is removed and replaced. SenSura® 1-piece urostomy is available with flat or convex baseplates. The broad assortment consists of pouches in a variety of sizes available in transparent or opaque colour.



## **SenSura® 2-Piece Urostomy**

### **SenSura® Flex**

SenSura® Flex is a two-piece appliance with adhesive coupling, providing the benefits of a one-piece – the unique combination of security and flexibility. The pouch is fixed securely to the SenSura® double-layer adhesive baseplate by a strong adhesive coupling system. The adhesive coupling is flexible, allowing freedom of movement, making SenSura® Flex comfortable and discreet to wear. The pouch can be removed from the SenSura® double-layer baseplate, allowing the pouch to be changed more often than the adhesive baseplate.



### **SenSura® Click**

SenSura® Click urostomy has a range of features designed to offer comfort and discretion. A non-return valve prevents the urine from washing over the stoma, reducing the risk of urinary infections. It has a soft tuck-away outlet which is easy to use and can be tucked away when not in use. This prevents it from rubbing on your leg, offering more comfort and discretion. It is divided into multichambers which spread the urine evenly within the pouch. This allows it to have a flat and discreet profile and reduces splashing sounds.





# SenSura<sup>®</sup> Baseplates

## **SenSura<sup>®</sup> Standard Baseplate**

The original adhesive for optimal skin friendliness.

Available in SenSura<sup>®</sup> Click and Flex.



## **SenSura<sup>®</sup> Convex Baseplate**

The oval shape of the SenSura<sup>®</sup> Convex baseplate is designed to help a stoma that is difficult to manage, for example, a stoma that is flush, retracted, or located in a skin fold. It applies light and directed pressure on the peristomal skin to reduce the risk of leakage by allowing the stoma to protrude. Available in SenSura<sup>®</sup> Click and Flex.



















































## **SenSura<sup>®</sup> Xpro Baseplate**

SenSura<sup>®</sup> Xpro is more erosion resistant than SenSura<sup>®</sup>. It is composed of adhesive which provides a fast and secure adherence that stays intact when exposed to aggressive output. It also protects the skin from aggressive output. Available in SenSura<sup>®</sup> Flex.



# SenSura® Product Overview

Product Code	Product Description	Size	Colour	Pack Size
<b>SenSura® Pouches (One Piece)</b>				
<b>SenSura® Closed</b>				
154500	SenSura 1-pc midi Closed C/Cut	10-66 mm	Opaque	30
154800	SenSura 1-pc maxi Closed C/Cut	10-76 mm	Opaque	30
154700	SenSura 1-pc maxi Closed C/Cut	10-76 mm	Transparent	30
154830	SenSura 1-pc maxi Closed C/Cut	35-76 mm	Opaque	30
154840	SenSura 1-pc maxi Closed C/Cut	35-76 mm	Opaque	30
<b>SenSura® Drainable</b>				
155800	SenSura 1-pc maxi Open C/Cut	10-76 mm	Opaque	30
155700	SenSura 1-pc maxi Open C/Cut	10-76 mm	Transparent	30
155090	SenSura 1-pc maxi Open	15-43 mm	Opaque	10
152060	SenSura 1-pc maxi Open	15-43 mm	Transparent	10
<b>SenSura® URO</b>				
118010	SenSura URO 1-pc Flat Bag, Midi	10-66 mm	Opaque	20
118020	SenSura URO 1-pc Flat Bag, Maxi	10-76 mm	Opaque	20
118030	SenSura URO 1-pc Flat Bag, Midi	10-66 mm	Transparent	20
118040	SenSura URO 1-pc Flat Bag, Maxi	10-76 mm	Transparent	20
118140	SenSura URO 1-pc Convex Light Bag, Maxi	50 mm	Transparent	10
118150	SenSura URO 1-pc Convex Light Bag, Maxi	60 mm	Transparent	10
<b>SenSura® CLICK Baseplates (Two Piece)</b>				
<b>SenSura® CLICK Standard Baseplate</b>				
100110	SenSura standard wear base plate 40mm	10-35 mm	-	5 
100210	SenSura standard wear base plate 50mm	10-45 mm	-	5 
100310	SenSura standard wear base plate 60mm	10-55mm	-	5 
<b>SenSura® CLICK Convex Baseplate</b>				
110210	SenSura standard Convex light 50mm	15-33 mm	-	5 
110310	SenSura standard Convex light 60mm	15-43 mm	-	5 
<b>SenSura® CLICK Pouches (Two Piece)</b>				
<b>SenSura® CLICK Closed</b>				
101540	SenSura 2pc midi closed	40mm	Opaque	30 
101550	SenSura 2pc midi closed	50mm	Opaque	30 
101560	SenSura 2pc midi closed	60mm	Opaque	30 
101650	SenSura 2pc maxi closed	50mm	Opaque	30 
101660	SenSura 2pc maxi closed	60mm	Opaque	30 

Product Code	Product Description	Size	Colour	Pack Size
<b>SenSura® CLICK</b>				
<b>SenSura® CLICK Drainable</b>				
103540	SenSura 2-pc midi Open	40 mm	Opaque	30 
103550	SenSura 2-pc midi Open	50 mm	Opaque	30 
103560	SenSura 2-pc midi Open	60 mm	Opaque	30 
103640	SenSura 2-pc maxi Open	40 mm	Opaque	30 
103650	SenSura 2-pc maxi Open	50 mm	Opaque	30 
103660	SenSura 2-pc maxi Open	60 mm	Opaque	30 
103840	SenSura 2-pc maxi Open	40 mm	Transparent	30 
103850	SenSura 2-pc maxi Open	50 mm	Transparent	30 
103860	SenSura 2-pc maxi Open	60 mm	Transparent	30 
<b>SenSura® Click URO</b>				
118440	SenSura URO 2-pc maxi	40 mm	Opaque	30 
118450	SenSura URO 2-pc maxi	50 mm	Opaque	30 
118460	SenSura URO 2-pc maxi	60 mm	Opaque	30 
118510	SenSura URO 2-pc midi	40 mm	Transparent	30 
118520	SenSura URO 2-pc midi	50 mm	Transparent	30 
118540	SenSura URO 2-pc maxi	40 mm	Transparent	30 
118550	SenSura URO 2-pc maxi	50 mm	Transparent	30 
118560	SenSura URO 2-pc maxi	60 mm	Transparent	30 
<b>SenSura® FLEX Baseplates (Two Piece)</b>				
<b>SenSura® FLEX Standard and Xpro Baseplate</b>				
101010	SenSura Flex Baseplate	35 mm	-	10 
101020	SenSura Flex Baseplate	50 mm	-	10 
101030	SenSura Flex Baseplate	70 mm	-	10 
101080	SenSura Flex Baseplate, standard	90 mm	-	5 
101040	SenSura Flex Xpro Baseplate	35 mm	-	10 
101050	SenSura Flex Xpro Baseplate	50 mm	-	10 
101060	SenSura Flex Xpro Baseplate	70 mm	-	10 
<b>SenSura® FLEX Convex Baseplate</b>				
113010	SenSura Flex Convex Baseplate	40 mm	-	5 
113030	SenSura Flex Convex Baseplate	50 mm	-	5 
113070	SenSura Flex Convex Baseplate	60 mm	-	5 
<b>SenSura® FLEX Pouches (Two Piece)</b>				
<b>SenSura® FLEX Closed</b>				
109040	SenSura Flex Maxi Closed Opaque	35 mm	Opaque	30 
109150	SenSura Flex Maxi Closed Opaque	50 mm	Opaque	30 
109180	SenSura Flex Maxi Closed Opaque	70 mm	Opaque	30 
<b>SenSura® FLEX Drainable</b>				
103220	SenSura Flex Maxi Open	50 mm	Transparent	30 
103270	SenSura Flex Maxi Open	70 mm	Transparent	30 
103900	SenSura Flex XXL Open	90 mm	Transparent	30 
103250	SenSura Flex Maxi Open	50 mm	Opaque	30 
103280	SenSura Flex Maxi Open	70 mm	Opaque	30 
103910	SenSura Flex XXL Open	90 mm	Opaque	30 
<b>SenSura® FLEX Urostomy</b>				
118270	SenSura Flex Uro Maxi	35 mm	Transparent	30 
118280	SenSura Flex Uro Maxi	50 mm	Transparent	30 

# *Alterna<sup>®</sup> Free*

*Spiral security with  
more freedom*

*Alterna<sup>®</sup> Free*


# Alterna<sup>®</sup> Spiral adhesive

The spiral adhesive is unique to Alterna<sup>®</sup>.  
The spiral form helps to ensure that  
Alterna<sup>®</sup> is both secure and skin-friendly.

**NEW**  
Increased  
flexibility

Skin-friendly spiral, keeps  
the skin comfortable and  
healthy


Adheres securely

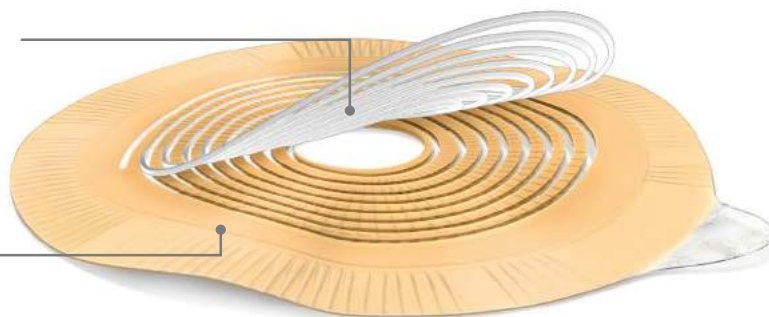

Now the spiral adhesive is enhanced to provide more comfort.

The adhesive now has extra  
flexibility for a more  
comfortable feel

Six flex-zones in the  
edge-area flex when the body  
does allowing for more  
comfortable movement

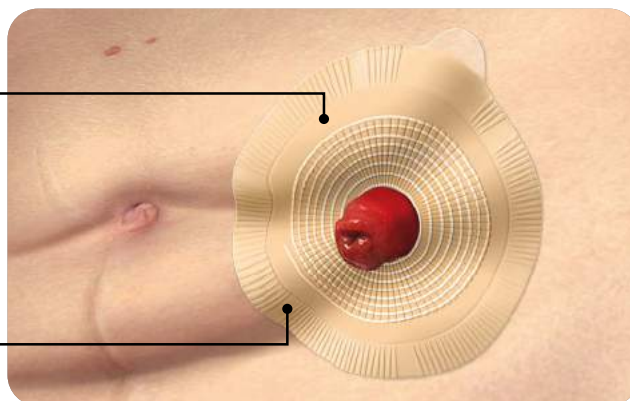



# Alterna<sup>®</sup> Closed Bags

## ***Alterna<sup>®</sup> 1-Piece Closed***

Alterna<sup>®</sup> 1-Piece Closed has an integrated multi-layer filter which helps avoid the risk of the pouch ballooning up and causing unnecessary bulges, and neutralises odour. The bag's soft backing fabric, in contact with the skin, features a lining made of soft material, allowing the skin to breathe and dries quickly after a shower or swim. The pouch is made of an odour-proof film.



## ***Alterna<sup>®</sup> 2-Piece Closed***

Alterna<sup>®</sup> 2-piece Closed offers a patented security ring with just a click to lock and unlock. When the system is locked, the pouch is safely and securely fastened to the barrier. The spiral-structure adhesive provides optimal adhesion and superior skin protection. The soft front makes the bag even more discreet. The bags feature an integrated filter, which allows air to escape. The filter is covered by a semi-permeable membrane, which allows air to pass, but not moisture. The membrane itself is covered by a net that allows air and moisture to pass, but not faeces. These two barriers around the active filter component provide double protection. The back of the bag, resting against the skin, features a lining made of extra soft material for comfort and added protection against skin irritation.



# Alterna® Free Drainable Bags

## **Alterna® Free 1-Piece Drainable**

Alterna® 1-Piece Drainable has a hideaway outlet, making the bag easy to empty and clean, and has an integrated Velcro® closure that can be easily tucked away – without the need for a closure clamp. An integrated three-layer filter avoids the risk of the pouch ballooning up and causing unnecessary bulges, and neutralises odour. The bag's soft backing fabric is strong and water repellent, allowing for easy drying after a wash or swim.



Also available without hideaway for use with a clamp.

## **Alterna® Free 2-Piece Drainable with Hideaway**

With Alterna® 2-Piece Drainable, the pouch is fixed securely to the baseplate by a mechanical coupling system. An audible 'click' confirms the secure locking of the pouch to the baseplate.

Before the pouch is locked into position, it can be rotated into the desired position, without having to remove the pouch from the baseplate. A self-adhesive closing clip is supplied with the pouch, giving a closure system that is both secure and easy to clean. The pouch can be removed from the spiral adhesive baseplate, allowing the pouch to be changed more often than the adhesive baseplate.



# Alterna® Urostomy Bags

## **Alterna® 1-Piece Urostomy Bag**

Alterna® 1-Piece Urostomy has a range of features designed to offer comfort and discretion. A secure and easy to use soft plastic outlet is easy to handle, and design of the pouch ensures there is no residual urine left in the pouch after emptying. The non-return valve prevents the urine from washing over the stoma, reducing the risk of urinary infections.



## **Alterna® 2-Piece Urostomy Bag**

Alterna® 2-piece Urostomy pouch integrates mechanical coupling with extra reassurance. The pouch is fixed securely to the baseplate by a mechanical coupling system. An audible 'click' confirms the secure locking of the pouch to the baseplate. The bag's secure and easy to use outlet is easy to handle, and design of the pouch ensures there is no residual urine left in the pouch after emptying. A non-return valve prevents the urine from washing over the stoma, reducing the risk of urinary infections.



# Alterna® Baseplates

## **Alterna® 1-Piece (Integrated Baseplate)**

Alterna® 1-piece appliances are available with flat or convex baseplates. The broad assortment consists of various hole sizes which can be customised, with pouches in a variety of sizes available in transparent or opaque colour.

Alterna® 1-piece is available with convexity, designed to help a stoma that is difficult to manage, for example a stoma that is retracted, or located in a deep skin fold or scar. It applies moderate to high pressure on the peristomal skin to reduce the risk of leakage by allowing the stoma to protrude. Ask your stoma care nurse for more information about convexity.

## **Alterna® 2-Piece Baseplates**

Alterna® 2-piece is available with a standard wear and extended wear baseplate, flat or with convexity.

**The broad range of baseplates are available in Custom-Cut:**

### **Custom-Cut**












































Custom-cut baseplates consist of a starter hole which can be customised and cut based on your own baseplate dimensions, ensuring a perfect and comfortable fit.



# Alterna® Product Overview

Product Code	Product Description	Size	Colour	Pack Size
<b>Alterna® One Piece Closed</b>				
174030	Alterna® 1-pc Closed bag, Maxi, with Filter (025160)	35 mm	Opaque	30
174040	Alterna® 1-pc Closed bag, Maxi, with Filter (025170)	40 mm	Opaque	30
174070	Alterna® 1-pc Closed bag, Maxi, with Filter	10-70 mm	Transparent	30
174040	Alterna® 1-pc Closed bag, Maxi, with Filter (025110)	20-55 mm	Opaque	30
174000	Alterna® 1-pc Closed bag, Maxi, with Filter (025110)	20-55 mm	Opaque	30
174020	Alterna® 1-pc Closed bag, Maxi, with Filter (025150)	30 mm	Opaque	30
174050	Alterna® 1- pc Closed bag, Maxi, with Filter	10-70 mm	Opaque	30
174100	Alterna® 1-pc Closed bag, Maxi, with Filter (024810)	20 -55 mm	Transparent	30
174200	Alterna® 1-pc Closed bag, Maxi, with Filter (024610)	20 -55 mm	Opaque	30
<b>Alterna® One Piece Drainable</b>				
174500	Alterna® 1-pc Open bag, Maxi, with Filter (025200)	10-70 mm	Opaque	30
174550	Alterna® 1-pc Open bag, Maxi, with Filter (024900)	10-70 mm	Transparent	30
174610	Alterna® 1-pc Open bag, Midi, with Filter (024500)	10-55 mm	Transparent	30
175000	Alterna® Free 1-pc Open bag with Hide-Away, Maxi, with Filter (138700)	12-70 mm	Opaque	30
175010	Alterna® Free 1-pc Open bag with Hide-Away, Maxi, with Filter (138600)	12-70 mm	Transparent	30
175110	Alterna® Free 1-pc Open bag, Convex, with Hide-Away, Maxi, with Filter (141060)	15-43 mm	Opaque	10
175150	Alterna® Free 1-pc Open bag, Convex, with Hide-Away, Maxi, with Filter (141030)	15-43 mm	Transparent	10
174700	Alterna® 1-pc Open bag, Convex, Maxi, with Filter (469480)	15-43 mm	Opaque	10
174710	Alterna® 1-pc Open bag, Convex, Maxi, with Filter (469380)	15-43 mm	Transparent	10
<b>Alterna® One Piece URO</b>				
174760	Alterna® 1-pc Urostomy bag, Maxi (055800)	10-55 mm	Opaque	30
174770	Alterna® 1-pc Urostomy bag, Maxi (055850)	10-55 mm	Transparent	30
174790	Alterna® 1-pc Urostomy bag, Midi (055750)	10-55 mm	Transparent	30
174810	Alterna® 1-pc Urostomy bag, Convex, Maxi (125950)	15-33 mm	Transparent	10
174820	Alterna® 1-pc Urostomy bag, Convex, Maxi (125960)	15-43 mm	Transparent	10
<b>Alterna® Paediatric</b>				
174670	Alterna® 1-pc Open bag, Children, with Filter	10-35 mm	Transparent	30
174800	Alterna® 1-pc Urostomy bag, Children (080090)	10-35 mm	Transparent	20



Product Code	Product Description	Size	Colour	Pack Size
<b>Alterna® Two Piece Base Plates</b>				
177000	Alterna® Base Plates for 2-pc bags (17700)	40 mm	-	5 
177050	Alterna® Base Plates for 2-pc bags (17750)	50 mm	-	5 
177090	Alterna® Base Plates for 2-pc bags (17800)	60 mm	-	5 
177160	Alterna® Base Plates for 2-pc bags (131710)	100/10-40 mm	-	5 
177190	Alterna® Base Plates for 2-pc bags (131810)	100/10-50 mm	-	5 
177200	Alterna® Base Plates for 2-pc bags (131910)	100/10-60 mm	-	5 
177020	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-40 mm	-	5 
177030	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-40 mm	-	5 
177070	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-50 mm	-	5 
177080	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-50 mm	-	5 
177110	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-60 mm	-	5 
177120	Alterna® Base Plates for 2-pc bags (Pre-Cut)	100/10-60 mm	-	5 
<b>Alterna® Two Piece Convex Base Plates</b>				
177360	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	40 mm (21)	-	5 
177380	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	50 mm (25)	-	5 
177400	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	50 mm (31)	-	5 
177420	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	60 mm (35)	-	5 
177460	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	50mm (15-33mm)	-	4 
177500	Alterna® Base Plates for 2-pc bags, Convex (Pre-Cut)	60mm (15-43mm)	-	4 
<b>Alterna® Two Piece Closed</b>				
017420	Alterna® 2-pc Closed, Mini Cap, with Filter	50 mm	-	30 
176000	Alterna® 2-pc Closed, Maxi, with Filter	40 mm	Opaque	30 
176010	Alterna® 2-pc Closed, Maxi, with Filter (016820)	50 mm	Opaque	30 
176020	Alterna® 2-pc Closed, Maxi, with Filter (016830)	60 mm	Opaque	30 
176030	Alterna® 2-pc Closed bag, Maxi, with Filter	40 mm	Transparent	30 
176040	Alterna® 2-pc Closed bag, Maxi, with Filter (016870)	50 mm	Transparent	30 
176050	Alterna® 2-pc Closed bag, Maxi, with Filter (016880)	60 mm	Transparent	30 
176060	Alterna® 2-pc Closed bag, Midi, with Filter (016710)	40 mm	Opaque	30 
176070	Alterna® 2-pc Closed bag, Midi, with Filter (016720)	50 mm	Opaque	30 
176080	Alterna® 2-pc Closed bag, Midi, with Filter (016730)	60 mm	Opaque	30 
<b>Alterna® Paediatric</b>				
139540	Alterna® Free 2-pc Open bag, Midi, with Filter	40 mm	Opaque	30 
139550	Alterna® Free 2-pc Open bag, Midi, with Filter	50 mm	Opaque	30 
139560	Alterna® Free 2-pc Open bag, Midi, with Filter	60 mm	Opaque	30 
139740	Alterna® Free 2-pc Open bag, Maxi, with Filter	40 mm	Transparent	30 
139750	Alterna® Free 2-pc Open bag, Maxi, with Filter	50 mm	Transparent	30 
139760	Alterna® Free 2-pc Open bag, Maxi, with Filter	60 mm	Transparent	30 
139840	Alterna® Free 2-pc Open bag, Maxi, with Filter	40 mm	Opaque	30 
139850	Alterna® Free 2-pc Open bag, Maxi, with Filter	50 mm	Opaque	30 
139860	Alterna® Free 2-pc Open bag, Maxi, with Filter	60 mm	Opaque	30 
<b>Alterna® Two Piece URO</b>				
176340	Alterna® 2-pc Urostomy bag, Maxi	40 mm	Opaque	20 
176350	Alterna® 2-pc Urostomy bag, Maxi	50 mm	Opaque	20 
176360	Alterna® 2-pc Urostomy bag, Maxi	60 mm	Opaque	20 
176370	Alterna® 2-pc Urostomy bag, Maxi (017570)	40 mm	Transparent	20 
176380	Alterna® 2-pc Urostomy bag, Maxi (017580)	50 mm	Transparent	20 
176390	Alterna® 2-pc Urostomy bag, Maxi (017590)	60 mm	Transparent	20 

# Hospital Assortment

SenSura® Mio

SenSura®

Alterna®







NEW

SenSura® Mio

# ***1-piece High Output appliance*** – for ileostomy patients and post-operative use

Only SenSura® Mio appliances feature BodyFit Technology® for **security** and **comfort**.

Now your patients can benefit from day one.

**SenSura® Mio**

# A broader range means everyone can find a SenSura® Mio that's *right for their needs*

Like all SenSura® Mio appliances, the new additions feature BodyFit Technology® to provide a secure and comfortable fit, no matter your patients' body shapes. The unique elastic adhesive follows the body's movements when bending and stretching, giving wearers the confidence to do more.



## The advantages of the NEW SenSura® Mio 1-piece High Output appliance:



Effective drainage through soft outlet. Connects to different night bags ensuring fewer appliance changes during the night.



Choose a flat or Convex baseplate, depending on the level of support your patient's stoma requires.



Available in both 70mm and 100mm for large and complicated stomas.



Scale on transparent bags provides visual indication of output volume in the bag.

## Why all SenSura® Mio appliances give patients the *confidence to do more*



### A unique elastic adhesive

The elastic adhesive follows the body's movements when bending and stretching. The adhesive also provides the tack, adhesion and water absorption to cope with more aggressive output.



### Discreet material

The smooth all-textile material makes it seem less like a medical appliance and more like a clothing item or accessory. Its specially chosen neutral colour is designed to stay unnoticed under all clothing colours – even white.



### A filter your patients can trust

Each appliance has a specially developed full-circle pre-filter that provides greater protection against clogging. It is proven to reduce incidents of ballooning by 61%.



NEW

SenSura® Mio

# ***2-piece High Output appliance*** – for ileostomy patients and post-operative use

Only SenSura® Mio appliances feature BodyFit Technology® for **security** and **comfort**.

Now your patients can benefit from day one.

**SenSura® Mio**



# A broader range means everyone can find a SenSura® Mio that's *right for their needs*

Like all SenSura® Mio appliances, the new additions feature BodyFit Technology® to provide a secure and comfortable fit, no matter your patients' body shapes. The unique elastic adhesive follows the body's movements when bending and stretching, giving wearers the confidence to do more.



## The advantages of the NEW SenSura® Mio 2-piece High Output appliance:




Available with a mini High Output bag for extra discretion when transitioning to the community.



Effective drainage through either drainable outlet or soft drainable outlet.



Choose a flat or Convex baseplate, depending on the level of support your patient's stoma requires.

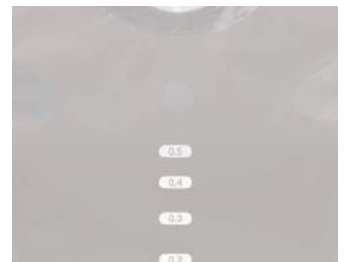

Scale on transparent bags provides visual indication of output volume in the bag.

## Click coupling



### Click – mechanical coupling

The Click coupling features a special lock-ring that makes a reassuring 'click' sound when correctly attached to the baseplate.

# SenSura® Post-Op bags

## ***Post-operative Care***

SenSura® post-op makes it possible for ostomates to feel secure right from the start, as it combines the SenSura® double-layer adhesive with the essentials in post-operative care.

## ***For easy access to the stoma***

Easy access is provided to the stoma through a window.

## ***For effective drainage***

SenSura® post-op bags allow for effective drainage through a soft outlet that connects easily to various night bags.



**SenSura® post-op is available in a variety of options to meet diverse post-operative needs:**

*Sterile and non-sterile pouches, and even a drainable pouch without a window.*

# Product Overview

## SenSura®

Product Code	Product Description	Size	Colour	Pack Size
<b>SenSura® Post Op</b>				
190200	SenSura Post-Op non-sterile with window	70 mm	-	6
190210	SenSura Post-Op non-sterile with window	100 mm	-	6
190300	SenSura Post-Op non-sterile without window	70 mm	-	6
190100	SenSura Post-Op sterile with window	70 mm	-	5
190110	SenSura Post-Op sterile with window	100 mm	-	5

## Alterna®

Product Code	Product Description	Size	Colour	Pack Size
<b>Alterna® Post Op Bag</b>				
128000	SenSura Post-Op non-sterile with window	10-100mm	-	6
128020	SenSura Post-Op non-sterile with window	10-100mm	-	6
128050	SenSura Post-Op non-sterile without window	10-70mm	-	6
128080	SenSura Post-Op sterile with window	10-70mm	-	5
128100	SenSura Post-Op sterile with window	10-100mm	-	5
128180	SenSura Post-Op sterile without window	10-70mm	-	5

## SenSura® Mio

Product Code	Product Description	Size	Colour	Pack Size
<b>SenSura® Mio One Piece High Output Post-Op Bag</b>				
186600	SenSura® Mio One Piece High-Output, Opaque with Inspection-Window POo	10-70mm	-	10
186610	SenSura® Mio One Piece High-Output, Opaque with Inspection-Window POo	10-100mm	-	10
186650	SenSura® Mio One Piece High-Output, POo Transparent	10-70mm	-	10
186660	SenSura® Mio One Piece High-Output, POo Transparent	10-100mm	-	10
				10
186700	SenSura® Mio One Piece High-Output, Soft Convex, Opaque with Inspection-Window POo	10-50mm	-	10
186710	SenSura® Mio One Piece High-Output, Soft Convex, POo Transparent	10-50mm	-	10
186750	SenSura® Mio One Piece High-Output, Light Convex, Opaque with Inspection-Window POo	10-43mm	-	10
186760	SenSura® Mio One Piece High-Output, Light Convex, POo Transparent	10-43mm	-	10
<b>SenSura® Mio Two Piece High Output Post-Op Bag</b>				
186010	SenSura® Mio Two Piece, Click High-Output, Opaque with Inspection-Window	50 mm		10
186020	SenSura® Mio Two Piece, Click High-Output, Opaque with Inspection-Window	60 mm		10
186030	SenSura® Mio Two Piece, Click High-Output, Opaque with Inspection-Window	70 mm		10
186110	SenSura® Mio Two Piece, Click High-Output, Transparent	50 mm		10
186120	SenSura® Mio Two Piece, Click High-Output, Transparent	60 mm		10
186130	SenSura® Mio Two Piece, Click High-Output, Transparent	70 mm		10
186400	SenSura® Mio Two Piece, Click High-Output, Transparent POo	50 mm		10
186410	SenSura® Mio Two Piece, Click High-Output, Transparent POo	60 mm		10
186420	SenSura® Mio Two Piece, Click High-Output, Transparent POo	70 mm		10
<b>SenSura® Mio Post-Op Bag</b>				
186800	SenSura® Mio Post-op, Sterile POo Transparent	10-70mm		5
186810	SenSura® Mio Post-op, Sterile POo Transparent	10-100mm		5
186900	SenSura® Mio Post-op, Non-Sterile POo Transparent	10-70mm		5
186910	SenSura® Mio Post-op, Non-Sterile POo Transparent	10-100mm		5

# Brava® supporting products

For customising appliance fit and for issues such as leakage and skin irritation.

**Brava®**



Sting-free  
easy  
removal



# Brava® Adhesive Remover

Quick and easy removal of adhesives

## The adhesive remover that makes a difference

Brava® Adhesive Remover lessens the chance of trauma on the skin by letting you and your patients easily remove the baseplate, leftover residue and any adhesive – without compromising the adhesion of their next appliance.

Brava® Adhesive Remover is available as a spray and also as a convenient wipe.

## Sting-free

As Brava® Adhesive Remover does not contain alcohol, it does not sting upon application.

## Skin friendly

Brava® Adhesive Remover contains silicone, which is gentle to the skin. Silicone-based adhesive remover also makes the baseplate easier to remove.

## Dries in seconds

No inconvenient delays before your patients are ready to apply a new baseplate.



87%

Less peel-force...<sup>1</sup>



...and tests suggest that the barrier function of the skin is maintained.<sup>1</sup>

Product name	Item number	Unit in package	Volume
Brava® Adhesive Remover spray	120105	1	50ml
Brava® Adhesive Remover wipe	120115	30	–

<sup>1</sup> The effect of Adhesive Remover Spray on skin stripping of ostomy products, 30-06-2016. Data on file.



Protects  
peristomal  
skin



## Brava® Powder

Absorbs moisture to lessen the risk of skin irritation

### Prevent skin moisture from causing problems

As well as causing a great deal of inconvenience and discomfort, skin moisture can also counteract the dry surface necessary for the adhesive to work properly. Brava® Powder absorbs moisture, keeping the skin dry and helping to stop skin irritation.

Brava® Powder can be used just around the stoma or under the adhesive.

### Protects peristomal skin

Brava® Powder reduces skin irritation related to moist skin.

### Helpful and friendly to skin

Brava® Powder helps to prevent red and irritated skin.



Product name	Item number	Unit in package	Size	Volume
Brava Powder	19075	1		25 g

Free of  
parabens,  
perfume and  
alcohol



# Brava® Barrier Cream

Gentle aid for dry or irritated skin

## Easy-to-use and long-lasting skin protection

Brava® Barrier Cream has a moisturising effect that helps to protect against skin irritation issues from when skin comes into contact with output, or when cleaning or removing the appliance.

Brava® Barrier Cream creates a long-lasting shield for the skin that ensures the baseplate adheres to the body. It also creates a long-lasting barrier that shields the skin from further damage from stoma output or adhesives.

## No parabens, perfume or alcohol

Brava® Barrier Cream is free of these potential irritants, making it gentle as well as effective.

## Does not compromise tapes or adhesives

Brava® Barrier Cream is designed to fit into your patients' appliance change routine. It allows tapes and adhesives to adhere so it does not need to be removed before applying the bag.

## Effective moisturising

Because Brava® Barrier Cream moisturises, it helps to manage sore, dry or irritated skin and prevent further damage.

Product name	Item number	Unit in package	Volume
Brava Barrier Cream	120005	1	60 ml

Provides  
a feeling of  
security



## Brava® Belt

Get a better fitting appliance – and more peace of mind

### The secure way to a better fit

The Brava® Belt discreetly provides extra security by helping to keep the baseplate in place.

Made of soft and comfortable material, the Brava® Belt comes in both standard and XL sizes that are adjustable to fit almost any body shape.



### Simple to wear

The Brava® Belt fits with all SenSura® and Assura® baseplates with belt ears.

There is a special belt for SenSura® Mio baseplates.

### Hygienic and reusable

The Brava® Belt is designed to be used multiple times and can be washed at 60 °C or higher.

### Help for all body shapes

The Brava® Belt can be particularly useful for those looking for the security of a more stable appliance.

Product name	Item number	Unit in package	Size
Brava Belt	04215	1	Standard
Brava Belt	04220	1	XL
Mio Belt	00423	1	Standard
Mio Belt	00424	1	XL

# Brava® Protective Seal

Easy-to-handle defence against leakage

NEW

## Feel secure with Brava® Protective Seal

If you have problems with leakage or skin irritation, the solution may be to improve the fit of your appliance with Brava® Protective Seal. Use it to create a snug fit around the stoma, or over dips and creases in the skin to seal any gaps between stoma and baseplate.

## Designed to protect against leakage

Seals gaps between your stoma and baseplate

The seal is easy to shape, so it fits snugly around the stoma and can be applied over uneven skin. It will stay firmly in place yet can be removed easily. Gives you dual-protection

The seal is output resistant so it will protect against the aggressive constituents of output. It also absorbs moisture, which helps to maintain skin health.

Try it and see what it can do for you!



Product name	Product Code	Pack Size	Thickness
Brava® Protective Seal, 18 mm	120350	30	2,5 mm
	120450	30	4,2 mm
Brava® Protective Seal, 27 mm	120370	30	2,5 mm
	120470	30	4,2 mm
Brava® Protective Seal, 34 mm	120390	30	2,5 mm
	120490	30	4,2 mm

# Brava® Mouldable Ring

Flexible and durable

## Less risk of leakage with a better fit

The most common issues that lead to leakage can be avoided with the Brava® Mouldable Ring. It gives a durable seal between the stoma and the baseplate, preventing output from undermining the adhesive.

As well as maintaining a better fit between the stoma and the baseplate, its durability protects the skin against stoma output.

The Brava® Mouldable Ring comes in thicknesses of 2.0 mm and 4.2 mm.

## Closes the gap for longer

Your patients can easily mould the Brava® Mouldable Ring into the required body shape so it fits snugly around the stoma and creates a seal between the stoma and baseplate. And because it's so durable it'll work for longer. Hinders loosening

The Brava® Mouldable Ring works excellently over dips, creases and folds – which helps create the even skin surface necessary to prevent the appliance from detaching.

## Protects and shields

The Brava® Mouldable Ring absorbs skin moisture and shields against output to keep your patients' skin healthy.



Product name	Item number	Unit in package	Size
Brava Mouldable Ring	120305	30	2 mm
Brava Mouldable Ring	120425	30	4.2 mm

Sting-free  
paste



# Brava® Paste

Kind on skin, tough on leakage

## Durable protection against leakage

Even if your patients have difficult folds and creases in their skin, Brava® Paste helps them to stay on top of leakage issues. Brava® Paste fills out difficult skin contours, maintaining a tight fit between the stoma and the appliance.

Brava® Paste is alcohol free so there is no itching or burning sensation when applying.

## No sting

Brava® Paste does not contain alcohol so there is no stinging sensation when it is applied.

## Shields skin

Brava® Paste helps to prevent skin issues by protecting against output. It also has a skin-friendly pH level.

## Resistant to breaking down or dissolving

Because Brava® Paste is so durable, it keeps working to help skin stay healthy.

Product name	Item number	Unit in package	Size	Volume
Brava Paste	120500	1		60 g

Easy  
to use



## Brava® Strip Paste

The filler that improves peristomal skin

### A tight seal despite scars and folds

The human body rarely provides stoma patients with the ideal surface for an appliance. Natural contours of the skin can lead to adhesion and leakage issues. Brava® Strip Paste can fill out even creases and deep folds, absorbing moisture and creating an even surface for the baseplate. This gives a better seal between the baseplate and skin, which helps reduce the risk of leakage.

### Flexible, durable and strong

Stoma users can be worried that their skin surface is too uneven or moist to create a strong seal. Brava® Strip Paste helps improve the peristomal skin. All that is necessary to do is to tear off the required length of strip and mould it into the required shape.

Product name	Item number	Unit in package
Brava Strip Paste	026555	10





# Brava® Elastic Tape

Appliance security for all shapes and sizes

## Soft elastic tape gives extra security by keeping the appliance more stable

Skin creases, folds and moisture all have a negative effect on keeping your patients' baseplate in place. By following body shape and movement, Brava® Elastic Tape keeps the baseplate secure. It also prevents the edges of the baseplate from lifting and rolling.

Brava® Elastic Tape comes in a variety of shapes to fit various body profiles and types of adhesives. The range also includes Brava® Elastic Tape Straight for square baseplates and Brava® Elastic Tape Belt, which provides extra support for outward areas such as a hernia.



## More freedom

Brava® Elastic Tape is developed with new BodyFit Technology®. This means greater security, longer wear time by preventing the baseplate edges from lifting – and more freedom of movement for your patients.

## Durable and skin friendly

Brava® Elastic Tape absorbs moisture, which makes the baseplate adhesive less prone to detach from the skin.

## Flexible usage

Patients can use one, two or more pieces depending on how much security and support they require.

Product name	Item number	Unit in package
Brava Elastic Tape	120700	20
Brava Elastic Tape Belt	120720	20
Brava Elastic Tape Straight	120740	20

Skin  
friendly



# Brava® Protective Sheet

An extra layer of skin protection

## Gentle peristomal skin protection

Patients with moist skin are often vulnerable to skin maceration and its attendant discomforts. The skin-friendly Brava® Protective Sheet is used under the baseplate, providing extra absorption in addition to the existing ostomy solution.

The Brava® Protective Sheet comes in a variety of sizes to fit different body shapes and needs.

## Creates a solid foundation

The Brava® Protective Sheet helps to avoid skin issues by keeping skin dry and providing a ready surface for applying a new appliance.

Product name	Item number	Unit in package	Size
Brava Protective Sheet	032105	10	10x10 cm
Brava Protective Sheet	032155	5	15x15 cm
Brava Protective Sheet	032205	5	20x20 cm

**Dries in  
seconds**



# Brava® Skin Barrier

## *An invisible friend for the skin*

### Shields against output and adhesives

Leakage and adhesives can cause patients a variety of skin problems. Brava® Skin Barrier leaves a thin breathable film on the skin that protects against potential damage but does not affect the adhesion of the next appliance.

The Brava® Skin Barrier is a silicone-based product so it does not sting when applied. It comes as a spray and a practical travel-friendly wipe.

### Prevents residue build-up

When your patients use Brava® Skin Barrier, residue either comes off when removing the adhesive or by gentle cleaning. This means it does not build up on the skin over time.

### Dries in seconds

With Brava® Skin Barrier, your patients will be ready to apply a new adhesive almost immediately.

Product name	Item number	Unit in package	Volume
Brava Skin Barrier spray	120205	1	50 ml
Brava Skin Barrier wipe	120215	30	-

Odour  
reduction

# Brava® Lubricating Deodorant

## Multi-lubricating effect

### Deodorant and lubrication in one

Odour from the output or output gathering around the stoma is a challenge for many patients. Specially designed to be used inside the bag, Brava® Lubricating Deodorant helps against odour when changing or emptying an appliance, while the multi-lubricating effect makes it easier for output to go to the bottom of the bag.

Brava® Lubricating Deodorant comes in a 240 ml bottle and an easy single-use sachet, which is ideal for patients who are away from home.



### Smooth removal and emptying

As well as making the bag easier to empty, Brava® Lubricating Deodorant's multi-lubricating effect also helps to keep the content at the bottom of the bag rather than around the stoma.

### Bag seal and adhesive still strong

Brava® Lubricating Deodorant is optimised for use with ostomy appliances so it does not degrade bag seals or affect the function of the adhesive and bag.

### Skin friendly

Although Brava® Lubricating Deodorant is designed to be applied inside the ostomy pouch, it is a skin-friendly product so there are no issues with it coming into contact with skin.

Product name	Item number	Unit in package	Size	Volume
Brava Lubricating Deodorant	120610	1		240 ml
Brava Lubricating Deodorant	120600	20	Small	7.5 ml



# Fistula and Wound Management System

## Environmental and disposal information

### Composition of product

- Ethylene vinylacetate - EVA
- Polyethylene - PE
- Styrene-ethylene-butadiene-styrene - SEBS
- Polyethylene foam - PE
- Polyvinylidene chloride - PVdC
- Polyester - PET
- Polyolefin - PO
- Adhesive
- Potato starch
- Pigment
- Activated carbon
- Methylparaben
- Ethylparaben

### Accessory - PUMP

- Polyoxymethylene - POM
- Polypropylene - PP
- Polyisoprene - PI

### Accessory - Cover sheet, tracing guide

- Ethylene vinylacetate - EVA
- Polyethylene - PE

### Accessory - CLAMP

- Polypropylene - PP

### Composition of packaging

The retail and transport boxes are made of cardboard.

- Carrugated cardboard
- Polyester - PET
- Polyethylene - PE

Product Code	Product Name	Pack Size
140100	Bed drainage bag	10
140150	Drain part	12
140200	Flexible lid, mini, with filter	10
140300	Flexible lid, midi, with filter	10
140400	Flexible lid, maxi, with filter	10
140550	Fistula and Wound Mng System, mini (with filter)	6
140650	Fistula and Wound Mng System, midi (with filter)	3
140750	Fistula and Wound Mng System, maxi (with filter)	3

# Notes

This image shows a full page of blank, lined paper. It features approximately 20 horizontal blue lines spaced evenly across the page, typical of notebook or primary writing paper. The lines are thin and light blue, set against a plain white background. There are no margins, text, or other markings on the page.







Coloplast® develops products and services that make life easier for people with very personal and private medical conditions.

Working closely with the people who use our products, we create solutions that are sensitive to their special needs. We call this intimate healthcare.

Our business includes ostomy care, urology and continence care and wound and skin care.

We operate globally and employ more than 8,500 people.

Coloplast SA, P.O. Box 10697, Fourways East, 2055, Tel: 0861 265 675, Fax/Orders: (011) 700-5003, Care number: 0861 61 2273 / 0861 61 CARE, Co. Reg. No. 94/08533/10.

[www.coloplast.co.za](http://www.coloplast.co.za) The Coloplast logo is a registered trademark of Coloplast A/S.  
© All rights reserved Coloplast A/S, 3050 Humlebæk, Denmark. ZA/09/2017-19/Corp/015

